



CPAP

Ascend[®] Nasal Mask

Ascend to a higher level of comfort and convenience

Experience a higher level of comfort and convenience with the Ascend[®] Nasal Mask, complete with our active headgear connector designed to minimize leaks for active sleepers.

Touchless Spacebar

Adding comfort with its contactless design, our headgear attaches to the touchless spacebar and provides additional stability by keeping the mask seal in place when the user moves.

Replaceable Cushion

Over time, regular wear and tear on the cushion of a respiratory mask can lead to discomfort and may reduce the effectiveness of the mask. Our replaceable cushions can be removed and cleaned, helping to maintain better hygiene and mask performance.

S M L

Quick-Release

Our quick release masks feature a slot on the side that allows the user to easily put on or take off the mask and headgear.

AIR[°]gel[®]


Sleepnet's AIR[°]gel[®] is an advanced cushion technology that conforms to the contours of each individual face creating a secure and comfortable seal.

360° Rotating Elbow

Designed to provide flexibility and comfort, the 360° rotating elbow allows the user to change sleep positions, move their head, or adjust the mask without discomfort.

Active Headgear Connector

Ideal for active sleepers, our active headgear connector is designed to move and adjust while maintaining a tight yet comfortable seal and preventing air leaks.

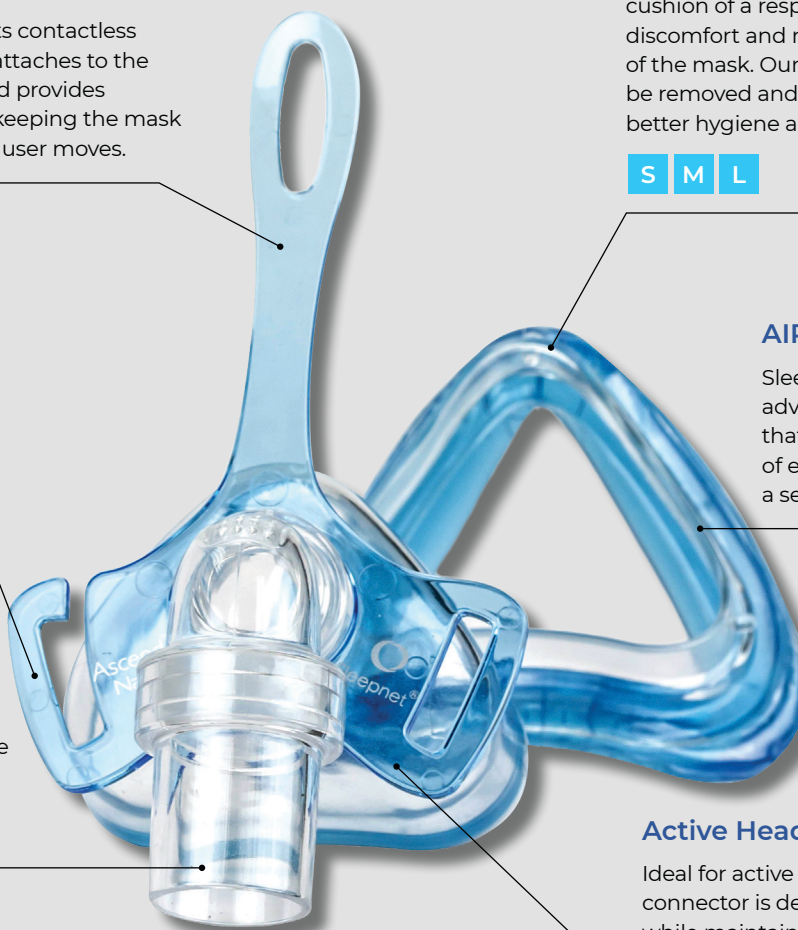



What makes our Ascend® Mask different?

Comfort - One of the most significant advantages of AIR°gel® technology is the unparalleled comfort it provides. The gel molds to the user's face, reducing pressure points and maximizing comfort.

Better Seal - the Ascend mask is designed to accommodate most facial features and provide superior stability while maintaining a comfortable seal.

Ease of Movement - Our active headgear connector is designed to move and adjust as the wearer moves, allowing the mask to maintain a seal and prevent air leaks. This can be especially helpful for active sleepers, allowing for ease of movement during sleep.

Available for resupply - Our replaceable cushions can be removed and cleaned, helping to maintain better hygiene and mask performance.

Easy-to-Replace Cushion S M L



Assembly: Place the top alignment pin in the notch. Push the bottom until all pins are secure.



Disassembly: Squeeze the sides of the cushion frame to release from mask.




Ascend® Nasal Mask with AFit Headgear Specifications:

Product Number	Description	Cushion Size	Mask Weight with Cushion
50174	Ascend® Nasal Mask system and headgear	S,M,L	N/A
50179	Ascend® Nasal Mask system with attached headgear	S,M,L	N/A
50157	Ascend® Nasal Mask with attached headgear	S	57.9 g
50158	Ascend® Nasal Mask with attached headgear	M	58.9 g
50159	Ascend® Nasal Mask with attached headgear	L	63.0 g
50187	Ascend® Nasal Mask without headgear attached	S	57.9 g
50188	Ascend® Nasal Mask without headgear attached	M	58.9 g
50189	Ascend® Nasal Mask without headgear attached	L	63.0 g
50201	AIR°gel cushion	S	NA
50202	AIR°gel cushion	M	NA
50203	AIR°gel cushion	L	NA
55058	AFit headgear only	N/A	N/A



Also by Sleepnet®
Ascend Full Face

Contact us today:

-  1-800-742-3646
-  info@sleepnetcorp.com
-  www.sleepnetmasks.com



The following product is protected by US patents 10,974,007, 9,199,054, and other patents pending



FDA
Cleared

Does not contain natural rubber latex

M90050 Version 0

