



CPAP & NIV

Mojo[®] 2 Full Face Mask

Personalized Comfort

The Mojo[®] 2 Full Face Mask is a lightweight mask with adjustable four-point headgear designed to provide a personalized, comfortable fit during patient therapy. It also accommodates for changes in weight and facial asymmetry.

Available Sizes:

S M L XL



Custom Fit Technology


Sleepnet's Custom Fit Technology is our proprietary technology made from a flexible, pliable material that responds to gentle manipulation and accommodates changes in weight and facial asymmetry.

Oxygen Port

Allows easy delivery of supplemental oxygen.

Available Elbow Options

- Vented
- Non-Vented
- Non-Vented AAV

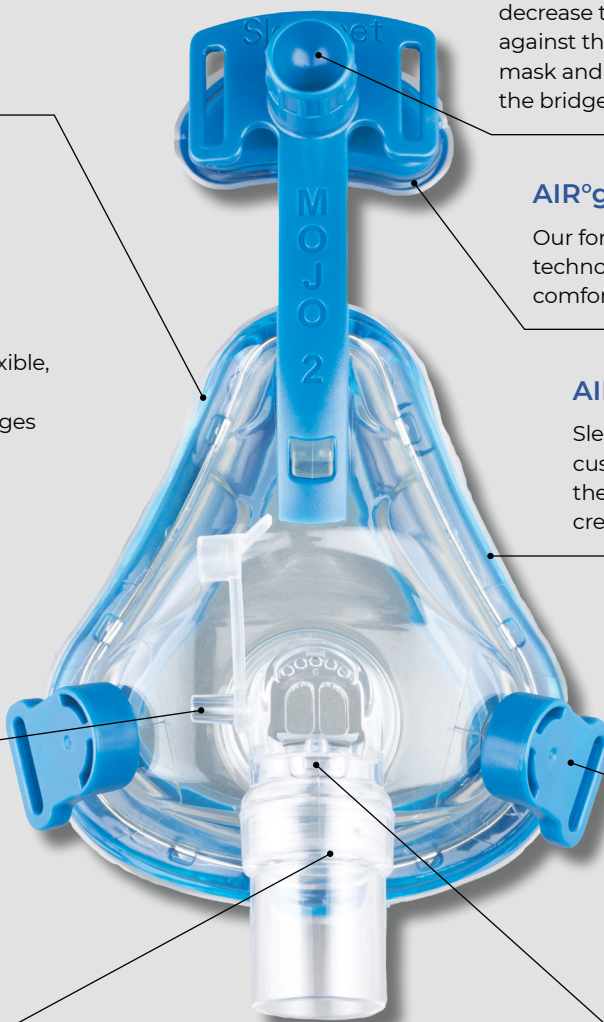

Adjustment Dial

By turning the dial, users can increase or decrease the pressure exerted by the pad against their forehead to help stabilize the mask and redistribute pressure away from the bridge of the nose.

AIR[°]gel[®] Forehead Cushion

Our forehead cushion uses AIR[°]gel[®] technology to provide enhanced comfort for the wearer.

AIR[°]gel[®]

Sleepnet's AIR[°]gel[®] is an advanced cushion technology that conforms to the contours of each individual face creating a secure and comfortable seal.

Easy-Release Magnetic Clips

Our magnetic clips are designed to simplify the process of attaching and detaching the mask from the headgear.

360[°] Rotating Elbow

Allows user to change sleep positions, move their head, or adjust the mask without discomfort.



What makes our Mojo[®] 2 different?

Comfort - Sleepnet's AIR^ogel[®] is a leading edge innovation developed and designed to enhance the comfort and effectiveness of our respiratory masks.

Personalized fit - Custom Fit Technology allows the user to easily shape the mask by gently bending and molding it along the sides, top, and bottom to provide a consistent and secure fit.

Adjustable - The Mojo[®] is designed with an adjustable dial and AIR^ogel[®] Forehead Cushion that are exclusive to the mask.

Easy on and off - our Easy-release magnetic clips are a user-friendly innovation designed to simplify the process of attaching and detaching the mask from the headgear.

Color coded elbow options for easy identification to meet your homecare and hospital needs.



Vented Mask with an anti-asphyxia valve for use with CPAP or bi-level devices.



AAV Non-Vented Mask with an anti-asphyxia valve for use with single limb circuits.



Non-Vented Mask for use with dual limb circuits.




Mojo[®] 2 Full Face Mask Specifications:

Product Number	Description	Size	Leak Port	Entrainment Valve	Mask Weight
50817	Mojo [®] 2 Full Face Mask Vented with Headgear	S	✓	✓	133.8 g
50818	Mojo [®] 2 Full Face Mask Vented with Headgear	M	✓	✓	141 g
50819	Mojo [®] 2 Full Face Mask Vented with Headgear	L	✓	✓	144.5 g
50820	Mojo [®] 2 Full Face Mask Vented with Headgear	XL	✓	✓	146.5 g
50912	Mojo [®] 2 Full Face Mask Non-Vented with Headgear	S	N/A	N/A	133.2 g
50913	Mojo [®] 2 Full Face Mask Non-Vented with Headgear	M	N/A	N/A	140.4 g
50914	Mojo [®] 2 Full Face Mask Non-Vented with Headgear	L	N/A	N/A	143.9 g
50915	Mojo [®] 2 Full Face Mask Non-Vented with Headgear	XL	N/A	N/A	145.9 g
50941	Mojo [®] 2 Full Face Mask AAV Non-Vented with Headgear	S	N/A	✓	133.7 g
50942	Mojo [®] 2 Full Face Mask AAV Non-Vented with Headgear	M	N/A	✓	140.9 g
50943	Mojo [®] 2 Full Face Mask AAV Non-Vented with Headgear	L	N/A	✓	144.4 g
50944	Mojo [®] 2 Full Face Mask AAV Non-Vented with Headgear	XL	N/A	✓	146.4 g
55051	EZFit Headgear only	One Size	N/A	N/A	N/A



Also by Sleepnet[®]
Ascend Full Face

Contact us today:

-  1-800-742-3646
-  info@sleepnetcorp.com
-  www.sleepnetmasks.com



The following product is protected by US patents 10,974,007, 7,938,117, D892,304.



FDA
Cleared

Does not contain natural rubber latex

M90083-01 Version 0



Sleepnet