



CPAP

V3 Full Face Mask

Providing extended comfort without compromise

The V3 Full Face Mask provides a restful home care experience and ensures a quick fit and comfort for the user.

AIR^ogel[®]

Sleepnet's AIR^ogel[®] is an advanced cushion technology that conforms to the contours of each individual face creating a secure and comfortable seal.

Adjustment Dial

By turning the dial, users can increase or decrease the pressure exerted by the pad against their forehead to help stabilize the mask and redistribute pressure away from the bridge of the nose.

AIR^ogel[®] Forehead Cushion

Our forehead cushion uses AIR^ogel[®] technology to provide enhanced comfort for the wearer.

Quick-Release

Our quick release masks feature a slot on the side that allows the user to easily put on or take off the mask and headgear.

360° Rotating Elbow


Allows user to change sleep positions, move their head, or adjust the mask without discomfort.

Oxygen Port

Allows easy delivery of supplemental oxygen.

Available Sizes:

S M L XL

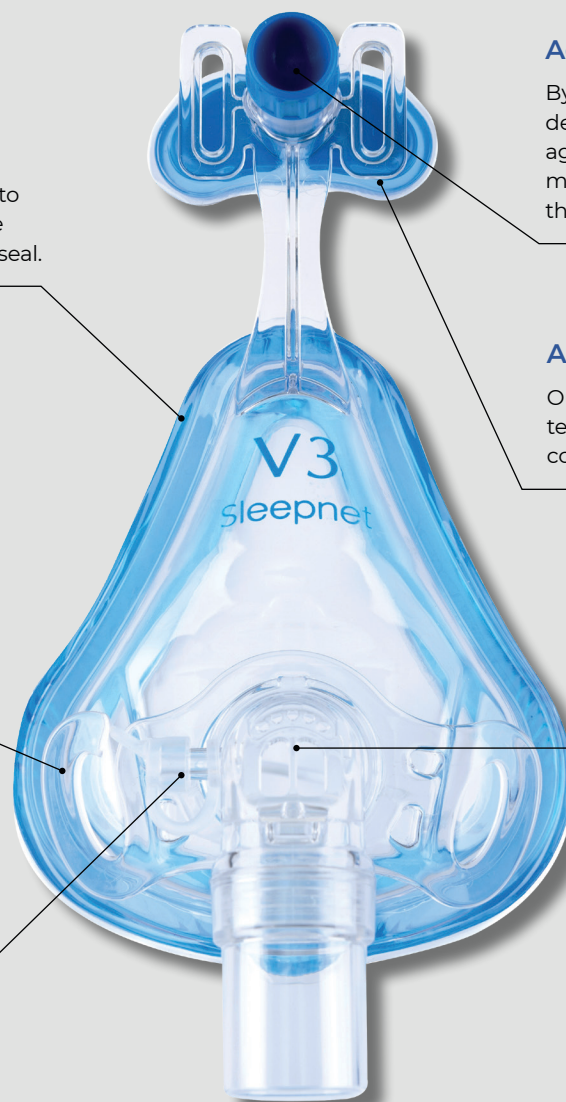



What makes our V3 different?

Comfort - Sleepnet's AIR°gel® is a leading edge innovation developed and designed to enhance the comfort and effectiveness of our respiratory masks.

Adjustable - The V3 is designed with an adjustable dial and AIR°gel® Forehead Cushion that are exclusive to the mask.

Less skin irritation - The V3 Full Face Mask is easily adjustable and built with Sleepnet's AIR°gel® cushion, which provides a comfortable seal while minimizing irritation.


Available in four sizes - Small, Medium, Large, and Extra Large.

The V3 Mask utilizes EZFit headgear for comfortable and secure treatment at night



V3 Full Face Mask Specifications:

Product Number	Description	Size	Leak Port	Entrainment Valve	Mask Weight
50863	V3 Full Face Mask Vented with Headgear	S	✓	✓	109.2 g
50864	V3 Full Face Mask Vented with Headgear	M	✓	✓	114.3 g
50865	V3 Full Face Mask Vented with Headgear	L	✓	✓	119.4 g
50866	V3 Full Face Mask Vented with Headgear	XL	✓	✓	120.1 g
55051	EZFit Headgear only	One Size	N/A	N/A	N/A

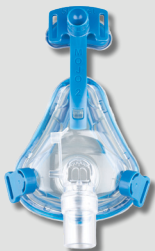

Also by Sleepnet®
Mojo® 2

Contact us today:

- ☎ 1-800-742-3646
- ✉ info@sleepnetcorp.com
- 🌐 www.sleepnetmasks.com



The following product is protected by US patents 10,974,007, 9,199,054, D892,304.



FDA
Cleared

Does not contain natural rubber latex

M90084-01 Version 0



Sleepnet